

MobiXRay

- Powerful
- Reliable
- Flexible
- Compact



MOBIXRAY_eng_Rev.2.1_02-17



Powerful, Efficient and Flexible



MobiXRay

Capacitor Assisted Mobile X-Ray



YOUR OEM PARTNER

SEDECAL
C/ Pelaya, 13 - Pol. Ind. "Río de Janeiro"
28110 - Algete (Madrid-SPAIN)
Central Offices:
Tel: (34) 91 - 628 0544
Fax: (34) 902 190 385
Email: info@sedecal.com

SEDECAL USA Inc
3190 N Konnicott Ave,
Arlington Heights, IL 60004
USA
Tel: 1 (847) 394 - 6960
Fax: 1 (847) 394 - 6966
Email: info@sedecalusa.com



MobiXRay


The mobile radiographic system **MOBIXRAY Capacitors** 16 kW, 20 kW or 32 kW, is an integrated capacitors Assisted RAD mobile. This system provides the solutions demanded by any Clinic and Hospital.

The system has the most innovative capacitor discharge High Frequency Generator that allows high output power from any standard power socket. Its excellent operability and maneuverability is ideal for any of the Hospital areas: surgical rooms, intensive care, emergency, etc.

Multiple positions for a quality performing

Max.
32kW

MobiXRay includes a powerful high frequency generator capacitor assisted, with a range of 16 kW, 20 kW, 32 kW Nominal power. The proven reliability of these generators assures constant performing during the whole day.

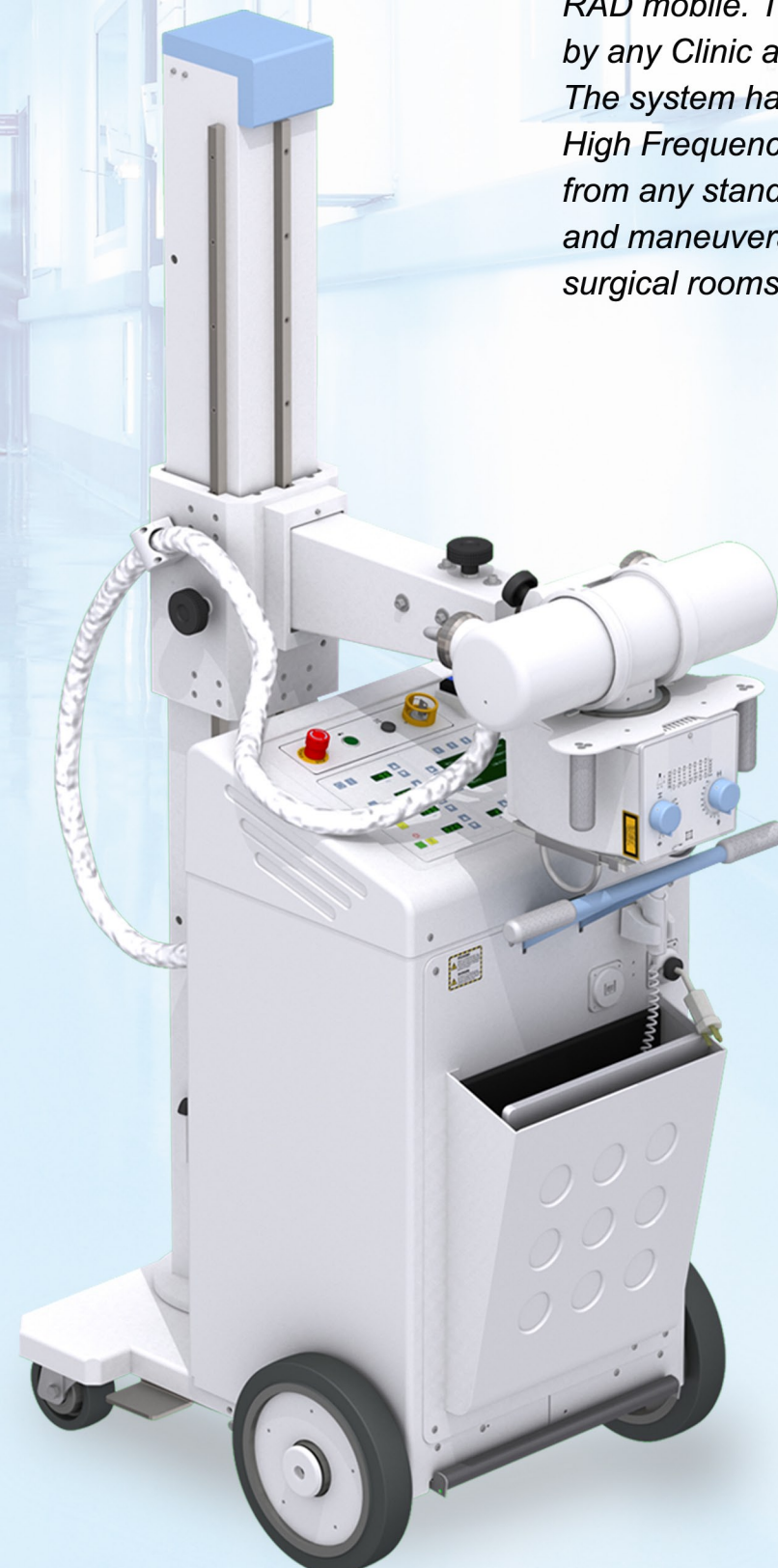


ADVANTAGES

TELESCOPIC ARM

The new system of telescopic arm has been designed with a smooth arm and maximum travel distance


Thanks to the **X-Ray tube rotation (180°)**, and the rotating collimator (90°), **MOBIXRAY** performs accurate exposures.

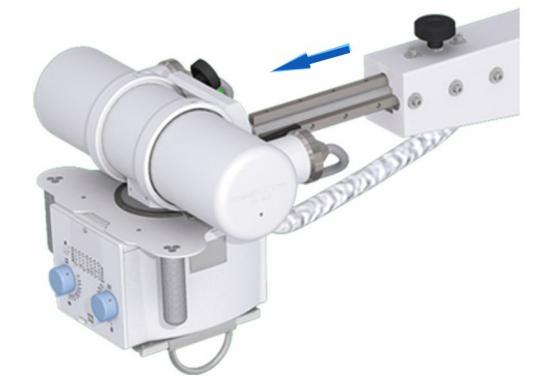

NEW ROTATING COLUMN

This column has been installed including an innovative system of rotation (conical bearing + axial bearing) that provides an accurate way of X-raying.

The **column rotates 315°**, which simplify x-raying for patients.

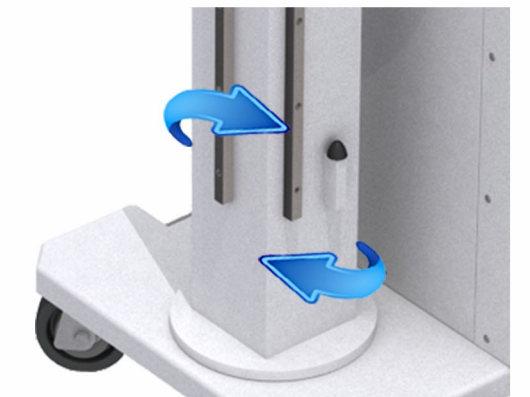

CONTROL CONSOLE & HANDLE

Large handle design to make driving easier, the design includes smooth hand grips, so handling the device is more comfortable.

The Control Console is very intuitive, also includes all the parameters for X-raying and the hand switch has been located at the front, to make the process easier.

